

## THE 2650KL

# DIGITAL WHEELCHAIR DUAL RAMP SCALE

## THE 2700KL

# DIGITAL WHEELCHAIR DUAL RAMP SCALE WITH EXTRA LARGE PLATFORM

The all-digital 2650KL and 2700KL support up to 1000 pounds and provide the optimum in quality, reliability, and ease of use. The ramps fold to save space, while the wide and low profile platform accommodates wheelchairs and provide a comfortable weighing experience for patients. The 2650KL features a medium-sized platform and folds for convenient storage. The 2700KL features the largest platform in the wheelchair collection to easily accommodate wide wheelchairs, a convenient choice for facilities where the scale remains in one location. Designed to meet the needs of hospitals, clinics, and acute and long-term care facilities, the scale's keypad features antimicrobial protection, a large non-skid low-profile platform, and dual ramps to easily wheel patients on and off the platform.

The remote display provides flexibility with the ability to be mounted on a wall or placed on a table. The scale's high-contrast color screen provides a user friendly interface that includes BMI, hold/release, reweigh, recall, zero, KG/LB, tare, pre-tare and more. The exclusive Everlock® function allows the measuring unit (LB or KG) to be permanently locked on the scale.

The scale provides easy integration and seamless communication with today's EMR applications, the 2650KL and 2700KL include a USB interface for connectivity with a computer, EMR software, or other electronic devices via various Health o meter® Professional Connectivity Solutions.



2700KL

## THERE'S STRENGTH IN OUR NUMBERS.

When you add up the numbers, it's not surprising that the Health o meter® Professional Scales is the industry-leading name in medical scales.

**120+** LOCAL  
FIELD SPECIALISTS

**100%**  
CONNECTIVITY

**#1** TRUSTED  
BRAND

**100** YEARS  
OF SERVICE

**Health o meter®**  
Professional

Health o meter® Professional 2650KL & 2700KL

weigh**easier**®

# TECHNICAL SPECIFICATIONS

## Capacity

1000 lb / 454 kg

## Resolution

0.2 lb / 0.1 kg

## EMR connectivity

Via USB, optional Pelstar® wireless technology

## Display

Digit Size:  $\frac{3}{8}$ " / 10 mm

Screen Size: 2  $\frac{3}{4}$ " x 1  $\frac{1}{2}$ " High-Contrast Color TFT-LCD Screen

## Power Source

100-240V adapter (ADPT30) included, 6 D batteries not included

## Functions

LB/KG conversion  
LB/KG lock out  
BMI  
Zero  
Tare  
Pre-tare  
Hold/release  
Reweigh  
Recall  
Auto zero  
Auto off  
Time/date  
Variable auto off time  
Audible/mute sound option  
Parent/child weighing function  
Everlock®

## 2650KL Platform Size (w x d x h)

31  $\frac{1}{2}$ " x 32  $\frac{1}{2}$ " x 2"  
(800 mm x 800 mm x 51 mm)

## 2700KL Platform Size (w x d x h)

43" x 42" x 2  $\frac{1}{2}$ "  
(1092 mm x 1067 mm x 65 mm)

## Ramp Size

2650KL: 10" x 31  $\frac{1}{4}$ " (254 mm x 794 mm)  
2700KL: 12  $\frac{3}{4}$ " x 41  $\frac{3}{4}$ " (324 mm x 1060 mm)

## 2650KL Product Footprint

51  $\frac{3}{4}$ " x 33  $\frac{5}{8}$ " x 5  $\frac{1}{4}$ "  
(1314 mm x 854 mm x 133 mm)

## 2700KL Product Footprint

69  $\frac{1}{8}$ " x 51  $\frac{1}{2}$ " x 3"  
(1756 mm x 1308 mm x 86 mm)

## Product Weight

2650KL: 82 lb / 37 kg  
2700KL: 187 lb / 85 kg

## Portability

2 wheels for ease of movement

## Warranty


2-year limited manufacturer warranty

# ORDERING INFORMATION

MODEL	DESCRIPTION
2650KL	Digital wheelchair dual ramp scale, LB & KG
2700KL	Digital wheelchair dual ramp scale with extra large platform, LB & KG
2650KG	Digital wheelchair dual ramp scale, KG only
2700KG	Digital wheelchair dual ramp scale with extra large platform, KG only
2650KL-BT	2650KL with built-in Pelstar® wireless technology
2700KL-BT	2700KL with built-in Pelstar® wireless technology
2650KG-BT	2650KG with built-in Pelstar® wireless technology, KG only
2700KG-BT	2700KG with built-in Pelstar® wireless technology, KG only
SS-2650KL	ScaleSurance® 2 year extended warranty for 2650 series*
SS-2700KL	ScaleSurance® 2 year extended warranty for 2700 series*
ADPT30	Adapter included, item # for replacement or spare
PLUGSET30	Plug adapters for use outside North America
BTPRINTERKIT-01	Wireless thermal printer kit

Scale also includes connectivity capabilities. To learn about connectivity options, visit [hom scales.com/ConnectivitySolutions](https://hom scales.com/ConnectivitySolutions)

To learn more, call 1-800-815-6615  
or visit [hom scales.com/2650KL](https://hom scales.com/2650KL)  
[hom scales.com/2700KL](https://hom scales.com/2700KL)

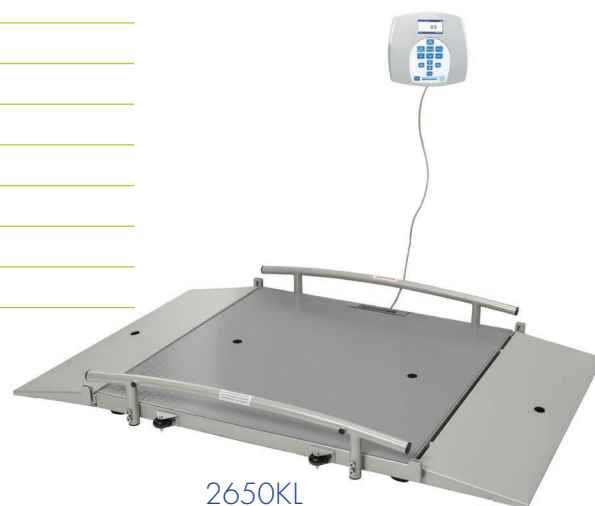

2650KL

**Health o meter®**  
Professional



\* Not available in all countries.