



Health o meter®
P r o f e s s i o n a l

weigh **easier**®
acute care
scales

Why Choose **Health o meter**[®] Professional Scales?

weigh**easier**[®]

What sets Health o meter[®] Professional apart is everything we do – from product design to technical support – is done with our customers in mind. Each decision and response is formed by first asking, “What will be weigh**easier** for our customer?” – whether the customer is a physician, healthcare system, distributor, or sales team.

Acute Care Scale Line

Dedicated to bringing you state-of-the-art products that are **Engineered with a Purpose**[®], the Health o meter[®] Professional hospital line was developed to improve workflow while keeping patients and caregivers safe. We are happy to support the market with best-in-class supply chain management and customer and technical support that has made us a leader in the marketplace.

#1 Medical Scale Brand in North America

Health o meter[®] has been a trusted brand since 1919. Today, Health o meter[®] Professional Scales is the U.S. market leader and the fastest-growing brand in Medical Scales.*

We offer:

- High-quality products that serve the entire continuum of medical care
- The industry’s largest network of trained and qualified field representatives
- Direct toll-free number for Technical Support with customer advocates based in the USA
- Fast shipments and industry-leading fill rates
- Standardization programs
- No-cost product evaluations

100% Focus on Healthcare

At Health o meter[®] Professional Scales, we design our products to help patients and the providers and institutions that care for them. We solely support the healthcare market and recognize the impact of healthcare reform, Value-Based purchasing, and other industry-specific initiatives. We understand providers are overwhelmed with patient care, focused on operating efficiently and monitoring costs. We are concerned with patient safety, supply chain efficiency, and technology standardization in response to our customer’s needs.

Products Designed in the USA

Health o meter[®] Professional’s US-based engineers and product specialists design high-value, high-quality devices for users around the globe.

Rated #1 by Medical Distributors & Healthcare Professionals!

Health o meter[®] Professional products are manufactured, designed and owned by Pelstar LLC. Health o meter[®] is a registered trademark of Sunbeam, Inc. used under license. ©2024 Pelstar[®], LLC

*As reported by Clarivate’s US Medical Supply Distribution data.

All Health o meter[®] Professional products are registered with the US Food & Drug Administration.

Page 4

**Heavy-Duty
Antimicrobial
Platform Scale**

Page 6

**Antimicrobial
High-Resolution
Digital Neonatal/
Pediatric Tray
Scale**



Page 8

**Antimicrobial
Wall-Mounted
Fold-Up
Wheelchair Scale**

Page 10

**Patient Transfer
Scale**



Page 11

**Digital Wet
Diaper/Lap
Sponge/Organ
Scale**

Page 12

**Automated
Calibrator**

Page 14

Connectivity

Page 15

Support

**ISO 9001:2015
ISO 13485:2016
CERTIFIED**

Health o meter®

P r o f e s s i o n a l



▶ 3105KL-AM



3105KL-AM

Heavy-Duty Antimicrobial Platform Scale

The 3105KL-AM has a robust, innovative design with enhanced safety and mobility features, including a low profile and large stable platform, an antimicrobial keypad and coating, smooth-rolling wheels, and an easy-grip handle. The scale makes it easy for clinicians to efficiently weigh patients and maintain their scale with multiple intuitive functions, a one-step calibration process, and easy to assemble instructions.

Capacity / Resolution

1000 lb (454 kg) / 0.2 lb (0.1 kg)

EMR connectivity

Via USB, optional Pelstar® wireless technology

Display

Digit Size: 3/8" / 10 mm
Screen Size: 2 3/4" x 1 1/2" High-Contrast
Color TFT-LCD Screen

Power Source

6 D batteries included, 100-240V adapter (ADPT30) not included



Platform Size (w x d x h)

19 3/4" x 21 1/4" x 2 1/8"
(502 mm x 540 mm x 54 mm)

Product Footprint (w x d x h)

23 1/4" x 28" x 52 1/4"
(591 mm x 711 mm x 1327 mm)

Product Weight

86 lb / 39 kg

Portability

2 smooth-rolling wheels and easy-grip handle

Warranty

2-year limited manufacturer warranty



| Model | Description |
|-----------------|--|
| 3105KL-AM | Antimicrobial heavy-duty platform scale, LB & KG |
| 3105KG-AM | Antimicrobial heavy-duty platform scale, KG only |
| 3105KL-AM-BT | 3105KL-AM with built-in Pelstar® wireless technology |
| 3105KG-AM-BT | 3105KG-AM with built-in Pelstar® wireless technology |
| 245EHR-3105 | Digital height rod for the 3105 series of scales |
| 201HR-3105 | Mechanical height rod for the 3105 series of scales |
| SS-3105KL-AM | ScaleSurance® 2 year extended warranty for 3105 series |
| ADPT30 | Optional adapter for 3105 series of scales |
| PLUGSET30 | Plug adapters for use outside North America |
| BTPRINTERKIT-01 | Wireless thermal printer kit |

Functions

- ✓ LB/KG conversion
- ✓ LB/KG lockout
- ✓ BMI
- ✓ Zero
- ✓ Tare
- ✓ Pre-tare
- ✓ Hold/release
- ✓ Reweigh
- ✓ Recall
- ✓ Auto zero
- ✓ Auto-off
- ✓ Time/date
- ✓ Variable auto-off time
- ✓ Audible/mute sound option
- ✓ Parent/child weighing function
- ✓ Everlock® - unit measure lock
- ✓ Enabling unit view with unit of measure locked

2210KL-AM / 2210KL4-AM

Antimicrobial High-Resolution Digital Neonatal/Pediatric Tray Scale

The 2210KL-AM neonatal scale was designed for use in the NICU, nursery, delivery room, and pediatric practices, where a high degree of accuracy is critical to the smallest of patients. The smooth-edged removable tray has a patent-pending mounting system that eliminates breakage, making it easy to clean. The scale will quickly lock despite the movement of the infant. The scale's backlit display makes reading a patient's weight easier at almost any angle and lighting condition. The high 1-gram resolution scale features a capacity of 20,000 g/20 kg/45 lb. The 2210KL-AM's user-friendly digital interface includes breast milk intake, auto-weight hold, hold/release, and more.

Capacity

45 lb / 20 kg / 20,000 g

Resolution

1 g (to 10,000 g);
5 g (10,001 to 20,000 g)

EMR connectivity

Via USB, optional Pelstar® wireless technology

Backlit LCD Display

5 ¼" x 1 ⅛" (133 mm x 28 mm)
4 levels of brightness

Power Source

6 D batteries and 100-240V adapter (ADPT30) included

Tray Size (w x d x h)

(2210KL-AM) 24" x 15 ¼" x 2 ¾"
(610 mm x 388 mm x 71 mm)
(2210KL4-AM) 26 ½" x 15 ¼" x 2 ¾"
(674 mm x 388 mm x 71 mm)

Product Footprint (w x d x h)

(2210KL-AM) 24 ¾" x 17 ¾" x 6"
(630 mm x 450 mm x 153 mm)
(2210KL4-AM) 26 ½" x 17 ¾" x 6"
(674 mm x 450 mm x 153 mm)

Product Weight

(2210KL-AM) 28 lb / 13 kg
(2210KL4-AM) 29 lb / 13.5 kg

Measuring Tape

0" - 24" / 0 cm - 60 cm

Portability

Easy movement between rooms is facilitated by durable handles and optional rolling carts

Warranty

2-year limited manufacturer warranty

| Model | Description |
|-----------------|--|
| 2210KL-AM | Antimicrobial high-resolution digital neonatal/pediatric tray scale, LB & KG |
| 2210KG-AM | Antimicrobial high-resolution digital neonatal/pediatric tray scale, KG only |
| 2210KL-AM-BT | 2210KL-AM with built-in Pelstar® wireless technology |
| 2210KG-AM-BT | 2210KG-AM with built-in Pelstar® wireless technology |
| 2210KL4-AM | Antimicrobial high resolution digital neonatal/pediatric 4-sided tray scale, LB & KG |
| 2210KG4-AM | Antimicrobial high resolution digital neonatal/pediatric 4-sided tray scale, KG only |
| 2210KL4-AM-BT | 2210KL4-AM with built-in Pelstar® wireless technology |
| 2210KG4-AM-BT | 2210KG4-AM with built-in Pelstar® wireless technology |
| 2210KL-AM-C | Antimicrobial high-resolution digital neonatal/pediatric tray scale, Canada only |
| 2210CART-AC | Multi-function Acute Care cart for 2210KL-AM/ 2210KL4-AM series of scales |
| 2210CART | Optional rolling cart for 2210 series |
| 2210EHR | Digital Height Rod for the 2210KL-AM Series of Scales* |
| SS-2210KL | ScaleSurance® 2-year extended warranty for 2210 series |
| ADPT30 | Adapter included, item # for replacement or spare |
| PLUGSET30 | Plug adapters for use outside North America |
| BTPRINTERKIT-01 | Wireless thermal printer kit |

Functions

- ✓ LB KG/G conversion
- ✓ LB/KG/G lockout
- ✓ Zero
- ✓ Tare
- ✓ Pre-tare
- ✓ Reweigh
- ✓ Recall
- ✓ Auto zero
- ✓ Auto-off
- ✓ Motion detection
- ✓ Everlock®
- ✓ Hold/release
- ✓ Variable auto-off time
- ✓ Audible/mute sound option
- ✓ Disable/enable backlight for display
- ✓ Enabling unit view with unit of measure locked
- ✓ Auto weight hold for 30 to 60 seconds
- ✓ Breast milk intake

*Available for use with 2210KL-AM series 2-sided tray only



▶ 2210KL4-AM



2900KL-AM

Antimicrobial Wall-Mounted Fold-Up Wheelchair Scale

The 2900KL-AM fold-up wheelchair scale has a capacity of 1000 pounds (454 kilograms), and a space-saving profile that only extends 3.5" from the wall. The large easy-to-lift platform accommodates most wheelchairs and provides a comfortable weighing experience. The dual ramps offer easy on, easy off access and the anti-slip caution tape delivers enhanced safety. Other safety features include antimicrobial protection, a raised edge, a recessed handle, and no pinch points.

Capacity / Resolution

1000 lb (454 kg) / 0.2 lb (0.1 kg)

EMR connectivity

Via USB, optional Pelstar® wireless technology

Display

Digit Size: 3/8" / 10 mm

Screen Size: 2 3/4" x 1 1/2" High-Contrast Color TFT-LCD Screen

Power Source

100-240V adapter (ADPT30) included, 6 D batteries not included



Platform Size (w x d x h)

29 1/2" x 32 1/2" x 2 3/4"
(749 mm x 825 mm x 70 mm)

Ramp Size (w x d)

6 1/4" x 34"
(159 mm x 864 mm)

Product Footprint Unfolded (w x d x h)

58" x 37" x 37"
(1473 mm x 940 mm x 940 mm)

Product Footprint Folded (w x d x h)

58" x 3 1/2" x 37"
(1473 mm x 89 mm x 940 mm)

Product Weight

91 lb / 41 kg

Warranty

2-year limited manufacturer warranty



Functions

- ✓ LB/KG conversion
- ✓ LB/KG lockout
- ✓ BMI
- ✓ Zero
- ✓ Tare
- ✓ Pre-tare
- ✓ Hold/release
- ✓ Reweigh
- ✓ Recall
- ✓ Auto zero
- ✓ Auto-off
- ✓ Time/date
- ✓ Variable auto-off time
- ✓ Audible/mute sound option
- ✓ Parent/child weighing function
- ✓ Everlock® - unit measure lock
- ✓ Enabling unit view with unit of measure locked

| Model | Description |
|-----------------|--|
| 2900KL-AM | Antimicrobial wall-mounted fold-up wheelchair scale, LB & KG |
| 2900KG-AM | Antimicrobial wall-mounted fold-up wheelchair scale, KG only |
| 2900KL-AM-BT | 2900KL-AM with built-in Pelstar® wireless technology |
| 2900KG-AM-BT | 2900KG-AM with built-in Pelstar® wireless technology |
| SS-2900KL-AM | ScaleSurance® 2-year extended warranty for 2900 series |
| ADPT30 | Adapter included, item # for replacement or spare |
| PLUGSET30 | Plug adapters for use outside North America |
| BTPRINTERKIT-01 | Wireless thermal printer kit |

PTS-1000KG Patient Transfer Scale

A transfer board with a built-in medical scale, the PTS-1000KG is designed to significantly reduce treatment times in stroke protocols or other time-sensitive medical conditions where precise weight measurement of immobile patients is crucial. The scale is lightweight, portable, and can be utilized as part of an existing patient transfer process without disrupting the patient flow process. Health o meter® Professional Scales warrants that the Patient Transfer Scale will be accurate to within 500 grams so long as the scale is properly used and maintained.

Capacity / Resolution

250 kg / 0.5 kg

Display

3" x 1" / 76 mm x 25 mm

Backlight LCD

Power Source

3.7V rechargeable battery

Weighing Mat Size (w x d)

71" x 21"

(1805 mm x 533 mm)

Product Footprint (w x d x h)

71" x 27 ½" x 1"

(1805 mm x 700 mm x 30 mm)

Product Weight

25 lb / 11.4 kg

Measuring Tape

Height guide shown on mat

100 cm – 180 cm shown at 5 cm

Warranty

2-year limited manufacturer warranty



Functions

- ✓ Zero
- ✓ Hold/release
- ✓ Auto off



Shorten Door-to-Needle Time Time Saves Lives

Always use a Slide Sheet when using the Patient Transfer Scale –



Check out the Single Patient Use and Reusable SlideMate™

| Model | Description |
|---------------|---|
| PTS-1000KG | Patient Transfer Scale, KG only |
| PTS-1000KL | Patient Transfer Scale, LB & KG |
| PTS-1000KL-C | Patient Transfer Scale, LB & KG, Canada only |
| PTS-1000KG-C | Patient Transfer Scale, KG only, Canada only |
| SS-PTS-1000KL | ScaleSurance® 2 year extended warranty for Patient Transfer Scale |

3400KL / 3401KL Digital Wet Diaper/Lap Sponge/Organ Scale

The 3400KL / 3401KL is a stable, stainless steel scale that has a removable weighing pan or tray that can be autoclaved for sterilization. It has a highly precise 0.05 oz / 1 g resolution and is easy to clean for sterile environments. Highlights include an easy-to-read backlit display, long-life rechargeable battery, and tare functionality.

Capacity / Resolution

176 oz (5000 g) / 0.05 oz (1 g)

Display

Digit Size: 5/8" (16 mm)

Screen Size: 3 7/8" x 1" (97 mm x 27 mm) LCD

Power Source

100 – 240V battery charger / 6 volt rechargeable battery (included)

Weighing Surface Dimensions- pan (w x d x h)

9 1/2" x 11 3/4" x 2 1/2"

(241 mm x 298 mm x 64 mm)

Weighing Surface Dimensions- tray (w x d x h)

7 3/4" x 9 3/4" x 5/8"

(197 mm x 248 mm x 16 mm)

Product Footprint- 3400KL (w x d x h)

10 7/8" x 13" x 8"

(276 mm x 330 mm x 203 mm)

Product Footprint- 3401KL (w x d x h)

10" x 11 3/4" x 5 3/8"

(254 mm x 298 mm x 137 mm)

Product Weight 3400KL

11 lb / 5 kg

Product Weight 3401KL

10 lb / 4.5 kg

Warranty

2-year limited manufacturer warranty

Functions

- ✓ Zero
- ✓ Tare
- ✓ Auto Zero
- ✓ Auto Off
- ✓ g/oz Conversion
- ✓ Low Battery Indicator
- ✓ Backlit Display



| Model | Description |
|-----------|--|
| 3400KL | Digital Wet Diaper / Lap Sponge / Organ Scale with Pan |
| 3401KL | Digital Wet Diaper / Lap Sponge / Organ Scale with Tray |
| 3400KG | Digital Wet Diaper / Lap Sponge / Organ Scale with Pan, KG only |
| 3400KG-C | Digital Wet Diaper / Lap Sponge / Organ Scale with Pan, KG only, Canada only |
| SS-3400KL | ScaleSurance® 2 year extended warranty 3400KL |
| SS-3401KL | ScaleSurance® 2 year extended warranty 3401KL |
| 3400PAN | 3 pack of stainless steel pans for 3400KL & 3401KL |
| 3400TRAY | 3 pack of stainless steel trays for 3400KL & 3401KL |
| 3400ADPT | Battery charger for 3400KL & 3401KL |



CALIBOT-AC/ CALIBOT-PC Automated Calibrator

The Automated Calibrator is a robust innovative calibration system that allows biomedical engineers and other medical technicians to calibrate medical patient scales quickly and safely without having to transport and stack physical calibration weights on the scale. Eliminating the cumbersome physical calibration weights from the calibration process reduces the time needed for routine scale maintenance and eliminates the safety risks associated with handling heavy dead weights.

Capacity / Resolution

500 lb (220 kg) / 0.01 lb (0.005 kg)

Platform Size CALIBOT-AC (w x d)

32" x 27"
(813 mm x 686 mm)

Platform Size CALIBOT-PC (w x d)

17" x 17"
(432 mm x 432 mm)

Product Weight CALIBOT-AC

97 lb / 44 kg

Product Weight CALIBOT-PC

65 lb / 30 kg

| Model | Description |
|------------|-------------------------|
| CALIBOT-AC | Acute Care Calibrator |
| CALIBOT-PC | Primary Care Calibrator |



Connectivity

Health o meter® Professional Scales offers connectivity solutions that transmit weight, height, and BMI measurements to vital signs monitors, PCs and EMR systems. Automatic data transfer reduces transcription errors and is the fastest, easiest way to capture patient information. Device connectivity can improve a healthcare facility's workflow efficiency, patient safety and EMR compliance efforts.

We offer wireless connectivity solutions that are included with the scale when purchased as a "BT" model or available as after-market add-ons to scales already deployed in the field.

Our scales are designed to interface with any EMR, and our list of connectivity partners is always growing. We will work with any EMR company at the request of our customers to facilitate direct interfaces with their systems. Among the systems we are compatible with are:



CEIBA HEALTH

IOMT + VIRTUAL CLINICAL CARE

CEIBA

Health o meter® Professional scale customers can seamlessly integrate their height and weight patient data into their EMR systems with CEIBA's device-agnostic IoMT interoperability. CEIBA Healthcare's IoMT eConnect Box and Cloud platform allow customers to swiftly and securely integrate their existing and new Health o meter® Professional scales to 20+ EMR or third-party systems, including Epic, Cerner, and PointClickCare.

Customized Solutions

Our Product Management Department can work directly with an organization to develop a solution to integrate Health o meter® Professional scales with a facility's system. Please contact Scott Gottman at sgottman@homscales.com, or call 708-297-6582.


To learn more about our Connectivity Solutions, visit:

homscales.com/ConnectivitySolutions


Global Support

Health o meter® Professional Scales has over 170 skilled professionals based in the United States and Canada who manufacture, sell, service, and support our products. Our global headquarters in McCook, IL is home to our highly trained Customer Service and Technical Support teams, along with our engineering, logistics, and management team. On-site sales and support are provided in partnership with MedTech/MedCare, an industry-recognized leader in medical sales throughout the US and Canada.

Our American-engineered products are designed and manufactured to meet the requirements of the International Organization for Standardization (ISO 9001:2015 & ISO 13485:2016).

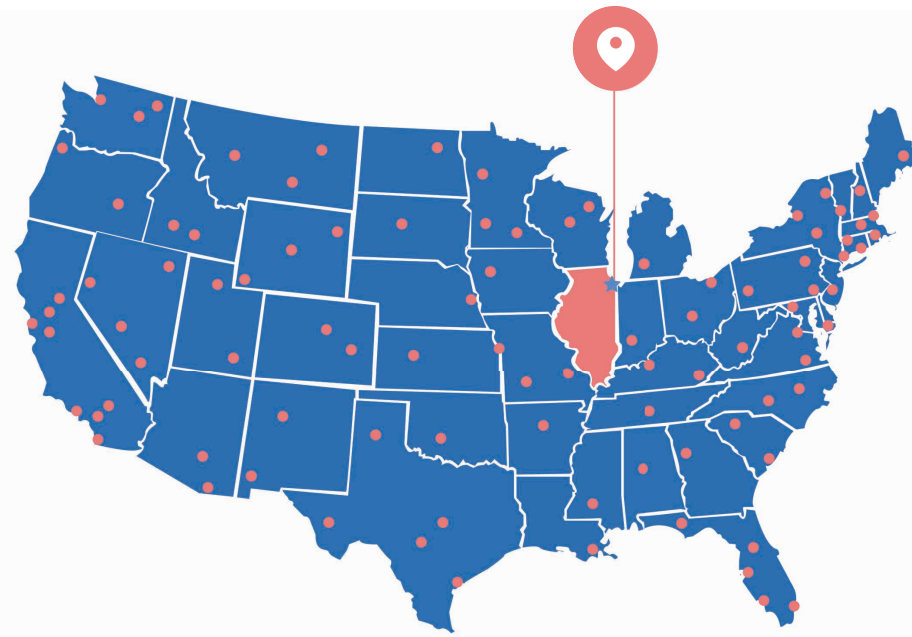


Health o meter® Professional Scales takes pride in designing world-class products for healthcare professionals around the globe. Our sales and support network includes select dealers and distributors throughout Latin America, Asia, Africa, and the Middle East.

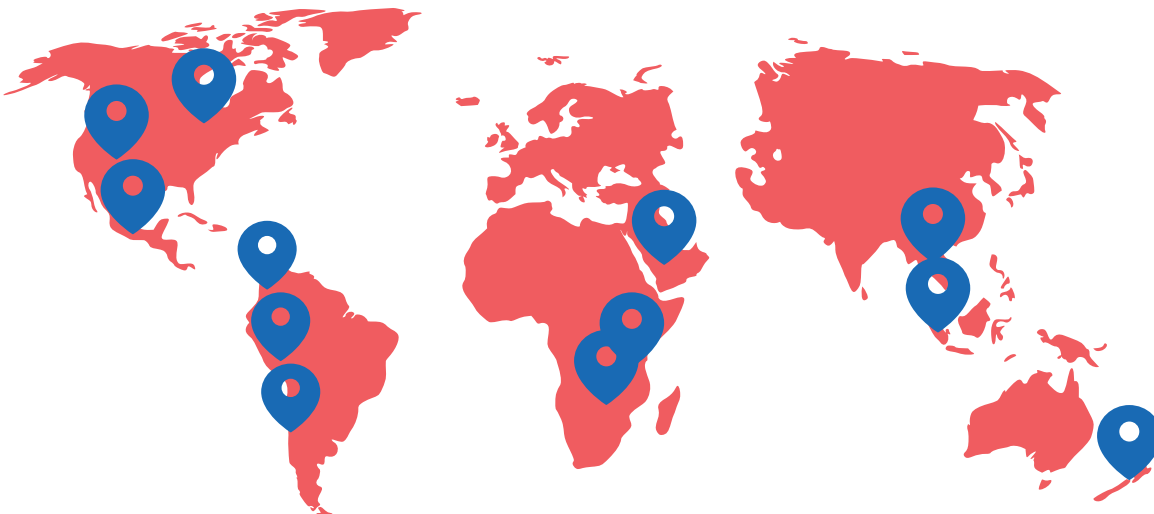

Health o meter[®]
P r o f e s s i o n a l

acute**care**

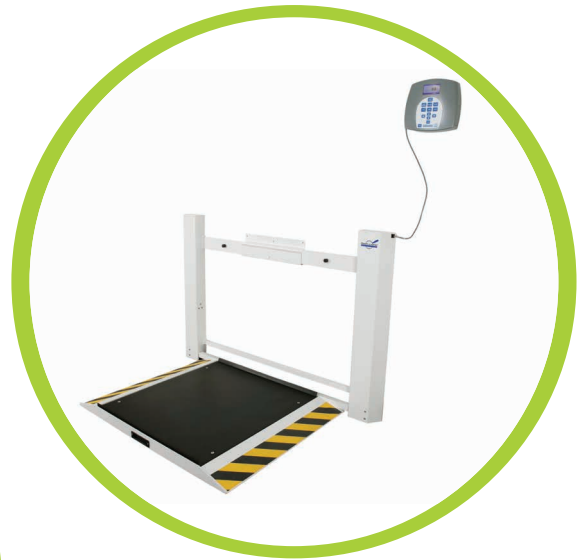

Contact Information

Main Number: 1 (708) 377-0600

Customer Service: 1 (800) 815-6615

Technical Support: 1 (800) 638-3722

Fax: 1 (708) 377-0601 or 1 (800) 654-7330

Email: homprocs@homscales.com

Corporate Headquarters:

Health o meter[®] Professional Scales / Pelstar, LLC

9500 West 55th St.

McCook, IL 60525-7110 USA